



Primipoesi



Dear customer,

Thank you for choosing this baby scale. Please read these instructions for use carefully and keep them for future use.

SPECIFICATIONS


- Accurately measures weight in 5g increments up to 20kg.
- HOLD feature captures precise weight measurement of an active baby.
- TARE feature subtracts added weight of diaper, clothing or blanket.
- Easy-to-read, white backlit LCD display.
- Unique touch button design.
- Includes 2xAAA batteries.
- Auto power off and low battery indication.

SCREEN



POWER SUPPLY

Please insert batteries or remove the insulating foil from the batteries before using. When "Lo" appears in the display, the batteries need to be replaced.

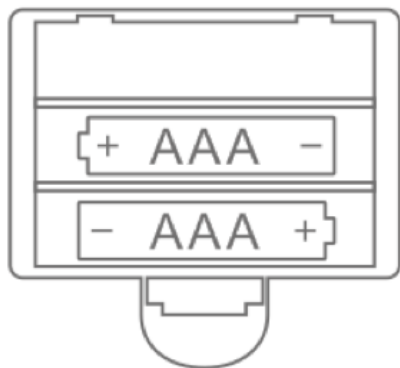

USAGE INSTRUCTIONS

The weight is displayed in steps of 5g from a minimum weight of 30g. Press the button (⏻) to switch on the scale.

When "0.000kg" appears on the display, the scale is ready for measurement. Place the baby on the scale carefully. Take particular care with the baby's head.

Hold Function:

The scale has a **HOLD** function: **AUTO HOLD**, **MECHANICAL HOLD**.

AUTO HOLD:

When the baby is placed on platform, they may move around. When they are still enough, the AUTO HOLD function will establish the weight on display. The sign (🔒) flashes and then weight is established.

MECHANICAL HOLD:

If the baby is moving around too vigorously, you may select the MECHANICAL HOLD function by pressing (🔒). The sign (🔒) flashes and weight is established on display. The sign (🔒) remains on display until you press (🔒) button or "TARE" button. Then it shows "0.000kg" which means the scale is ready for next weighing.

Tare Function:

If you prefer to place a towel or blanket on the scale and weigh the baby, please do so BEFORE switching the scale on. If the towel or blanket weighs over 30g, please press TARE button.

"TARE" is shown on display. Then place the baby onto towel or blanket. It shows baby's weight.

Units:

The scale is with units of kilograms (kg) and pounds (lb), which can be switched by pressing "UNIT" button.

SWITCH OFF

After weighing is finished, you can switch off the scale by long pressing (⏻) button. Or, the scale will automatically switch off if no weighing is carried out within 180 seconds.

CARE AND ATTENTION

- Place the scale on a hard, flat surface to ensure greatest accuracy.
- Cleaning: You can clean the scale with a damp cloth and a little washing up liquid if necessary. Never immerse the scale in water or rinse it under running water.
- Protect the scale against knocks, damp, dust, chemicals, drastic changes in temperature, electromagnetic field and source of heat that are too near.
- Keep packaging material away from children.
- DO NOT dispose of the old batteries in fire. The batteries may explode or leak. Recycle or manage them as hazardous waste.
- If the scale is not used in a long period, please take the batteries out.

WARRANTY - TERMS AND CONDITIONS

This product benefits of a warranty of 24 months on material and manufacturing defects, starting from the date of purchase (see sales slip).

The warranty of 24 months does not include damages caused by the usual use of parties identified as "consumable" (e.g., batteries, brush heads or parts subject to usury).

The legal guarantee of 24 months is void if:

1. The product has undergone aesthetic damage due to improper use not in accordance with instructions in the manual.
2. This product has been modified and/or tampered with.
3. The cause of the failure was due to poor maintenance of the individual components and/or accessories and/or supplies (e.g. oxidation and/or scaling due to the retention of water or other liquids, sediment blocking the sensor, leak of corrosive liquid from batteries).

The following is excluded from the legal guarantee of 24 months:

1. Costs related to replacements and/or repair of parts subject to wear or costs for ordinary maintenance of the product.
2. The costs and risks involved in transporting the product to and from the store where you purchased or otherwise authorized collection centre to receive the products under warranty.

OVERLOAD INDICATION

The scale can accept a maximum load of 20 kg. When overloaded, "EEEE" sign will appear on display.

3. Damage caused by or resulting from improper installation or improper use not in accordance with the directions in the instruction manual.

4. Damage due to natural disasters, accidental events or adverse conditions not compatible with the product.

5. Defects that have a negligible effect on product performance.

The manufacturer, distributor and all the parties involved in the sale do not assume any liability for losses and economic damage from any malfunction of the product. In accordance with current regulations the manufacturer, distributor and all the parties involved in the sale are not responding in any case for damages, including direct, indirect ones, loss of net income, loss of savings and additional damage and other details consequences going beyond the damage caused by the breach of warranty, contract, strict liability, wrongdoing or due to other causes, resulting from the use or inability to use the product and/or paper and electronic documents, including the lack of service.

For further information on the help service visit the website www.nuvitababy.com.