
1101

IT Set igiene orale
ES Kit de higiene bucal
EN Oral hygiene set
FR Kit d'hygiène bucco-dentaire
DE Mundhygiene-Set
PT Conjunto de higiene oral
GR Σετ στοματικής υγιεινής

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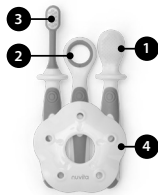EN

Thank you for purchasing the Nuvita 1101 Oral hygiene set. Please read the manual carefully before use and keep it for future reference.

IMPORTANT SAFETY INSTRUCTIONS

For the safety and health of you and your children, all the contents of this project are related to safety. Please be sure to abide by them.

WARNING: • This is not a toy. • Always use under adult supervision. • Do not allow children to use or chew on these items unsupervised. • The safety shield prevents the brush from going too far into the mouth. • Not suitable for boiling or hot steam sterilization.

INCLUDED IN THE SET: 1. Gum massager. 2. Tongue cleaner 3. Toothbrush. 4. Safety shield.

INSTRUCTIONS FOR USE:

GUM MASSAGER: • Always wash your hands before using. • Use to gently massage baby's gums and clean the tongue. • Use with water only, no toothpaste. **TONGUE CLEANER:** • Use to gently remove milk coating or food residues from the tongue. • Swipe from back to front. • Clean with warm water after use. **TOOTHBRUSH:** • Attach safety shield. • Apply a small amount of baby toothpaste. • Gently brush new teeth in circles. Encourage your child to brush with you. • Supervise at all times. **CLEANING, MAINTENANCE:** • Wash before and after each use with warm water and mild soap. • Do not boil or sterilize with hot water. • Air dry completely. • Store in a clean, dry place. **SAFETY GUIDELINES:** • Never leave the child alone with the brush. • Inspect all parts before each use. • Replace toothbrushes every 2-3 months or sooner if worn. •

Do not attach strings or cords. • Keep all items dry and clean. • Ensure all components are dry before placing in storage case. **WHEN TO REPLACE:** • At signs of wear, damage, or every few months as needed. • After illness, replace the brush to prevent reinfection.

WARNING! Clean gently with warm water and neutral detergent. Rinse well and dry thoroughly before storing in a cool, dry place. This product is not a toy. Always use under adult supervision. Discard at the first sign of deterioration, wear, or breakage.

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All rights reserved
Anteprima Brands International Ltd