
nuvita baby

SONIC CLEAN&CARE

- IT** Spazzolino elettrico per bambini da 3m+ fino a 5 anni
- ES** Cepillo de dientes eléctrico para bebés desde 3m+ hasta 5 años
- EN** Electric toothbrush for babies from 3m+ up to 5 years
- FR** Brosse à dents électrique pour bébés de 3m+ à 5 ans
- DE** Elektrische Zahnbürste für Babys von 3m+ bis 5 Jahre

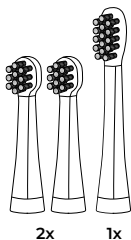- IT** Manuale di istruzioni
- ES** Manual de instrucciones de uso
- EN** Instruction manual
- FR** Manuel d'instructions d'utilisation
- DE** Bedienungsanleitung

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
Congratulations on your purchase, and welcome to Nuvita Baby! Discover all the compatible accessories and the Nuvita range of products on **www.nuvitababy.com**

The Nuvita Sonic Clean&Care Electric Toothbrush is specifically designed for babies and kids oral care. It has been approved by AIO, the Italian Dentist Association, and is in line with the World Health Organisation directives as regards to babies' oral care.

Regular tooth brushing is important for both babies and kids. It helps to remove the bacteria and plaque that cause tooth decay and gum disease. It is recommended by dentists that everybody brush their teeth twice a day – in the morning and before going to bed at night. Brushing also stimulates the gums, which helps to keep the teeth healthy and prevent gum disease. What starts out as a game with your little ones will turn into a healthy habit and help prevent cavities and oral care problems of a lifetime.

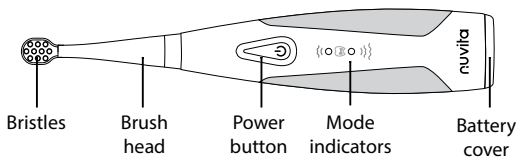
Content of the box:

1. Toothbrush body
2. 2 small and 1 medium brush heads
3. Instruction manual

Brushing Instructions:

1. Wet the bristles and apply appropriate amount of baby toothpaste if you wish (it can be used even without toothpaste)
2. Place the toothbrush bristles against the teeth at a slight angle towards the gum line.
3. Press the power button to turn on the toothbrush, gently move the brush. To select your desired speed and mode:
 - PRESS ONCE
GENTLE CLEANING
 - PRESS TWICE
DEEP CLEANING
 - PRESS THREE TIMES
PULSATING FUNCTION
4. Press the power button again to turn it off after finish brushing otherwise the baby toothbrush will turn off automatically in 2 minutes.
5. Last used speed/mode would be memorized for next time you power it on.
6. To avoid splashing, turn off your toothbrush before removing from the mouth by pressing the on / off button.

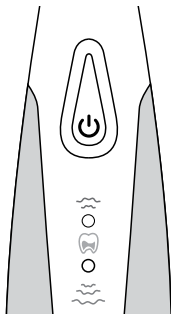Replacing the heads:

Dentists recommend replacing the brush heads maximum every 2 months to achieve optimal results. Use of damaged and ruined heads can damage gums and teeth. Use only the Nuvita Sonic Clean&Care replacement brush heads.

- Small size is suitable for age 3m+
(2 pcs included)
- Medium size is suitable for age 9m+
(1 pcs included)

Instructions:

1. Loosen and remove the installed head by turning it anti-clockwise then pulling it up (2).
2. Insert the new brush head into the notch and twist clockwise until the toothbrush head locks in place.

Cleaning:

Do not clean brush heads or handle in the dishwasher or a microwave.

Toothbrush handle


1. Remove the brush head and gently rinse the shaft area with warm water.
2. Use a damp cloth to wipe the entire surface of the handle.

Brush head

1. Rinse the brush head and bristles after each use.
2. Remove the brush head from the handle and rinse the brush head connection at least once a week with warm water.

Replacing the battery:

- Turn battery cover at the base of the toothbrush handle with a coin as direction (1).
- Insert a new AAA Alkaline battery with (+) facing outside.
- Install battery cover and turn it with a coin as (2).

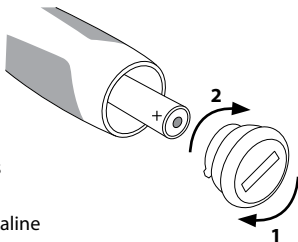Specifications:

Heads Bristles	
Material:	DuPont Nylon
Body:	ABS/TPE/PP
Power Supply:	1 X AAA battery
Vibration technology:	Sonic
Frequency:	16,000 / 22,000 strokes per minute selectable by user
Waterproof	
Battery Life:	up to 150 days
Product Size without brush head:	105 x 25 x 25 mm
Product Standard:	RoHS Directive 2011/65/EU 55014-1:2006+A1:2009 55014-2:1997+A1:2001+A2:2008

Environment:

Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment. Always remove the battery before discard. Dispose of the battery at an official collection point for batteries.

Warnings:

1. Always use under adult supervision.
2. Never leave the child alone with any of these accessories in his/her mouth or hand. This product is not a toy.
3. This product is designed for household use to clean teeth, gums and tongue only, do not use it for any other purposes.
4. Do not open, repair or alter the product without authorization.
5. Do not soak the product in water for long period of time.
6. Do not place in high temperature condition or directly under the sunlight.
7. Take out battery if no intention to use the product for an extended period of time.
8. Consult your dental professional before you use this product if your baby had oral or gum surgery
9. Use this product as described in this manual.
10. Do not use any other heads other than these manufactured under Nuvita brand and specifically made for this product.

WARRANTY - TERMS AND CONDITIONS

This product benefits of a warranty of 24 months on material and manufacturing defects, starting from the date of purchase (see sales slip). The warranty of 24 months does not include damages caused by the usual use of parties identified as "consumable" (e.g., batteries, brush heads or parts subject to usury). The legal guarantee of 24 months is void if: 1. The product has undergone aesthetic damage due to improper use not in accordance with instructions in the manual. 2. This product has been modified and/or tampered with. 3. The cause of the failure was due to poor maintenance of the individual components and/or accessories and/or supplies (e.g. oxidation and/or scaling due to the retention of water or other liquids, sediment blocking the sensor, leak of corrosive liquid from batteries). The following is excluded from the legal guarantee of 24 months: 1. Costs related to replacements and/or repair of parts subject to wear or costs for ordinary maintenance of the product. 2. The costs and risks involved in transporting the product to and from the store where you purchased or otherwise authorized collection centre to receive the products under warranty. 3. Damage caused by or resulting from improper installation or improper use not in accordance with the directions in the instruction manual. 4. Damage due to natural disasters, accidental events or adverse conditions not compatible with the product. 5. Defects that have a negligible effect on product performance. The manufacturer, distributor and all the parties involved in the sale do not assume any liability for losses and economic damage from any malfunction of the product. In accordance with current regulations the manufacturer, distributor and all the parties involved in the sale are not responding in any case for damages, including direct, indirect ones, loss of net income, loss of savings and additional damage and other details consequences going beyond the damage caused by the breach of warranty, contract, strict liability, wrongdoing or due to other causes, resulting from the use or inability to use the product and/or paper and electronic documents, including the lack of service. For further information on the help service visit the website

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All rights reserved
Anteprima Brands International Ltd