nuvita

MATERNO *SMART* *WEAR*

IT	Tiralatte elettrico indossabile
ES	Sacaleches eléctrico portátil
EN	Electric wearable breast pump
FR	Tire-lait électrique portable
DE	Elektrische tragbare Milchpumpe
PT	Tira-leite eléctrico vestível
GR	Φορητό ηλεκτρικό θήλαστρο
NO	Elektrisk bærbar brystpumpe
DA	Elektrisk bærbar brystpumpe

**Care
Flo®**
Technology

**Nip
Fit®**

1287CW

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
www.nuvitababy.com/warranty

ATTENTION

- Please do not leave this product in a high temperature environment for a long time.
- Please clean this product before and after use.
- Before breast pumping, please check again whether the related accessories are correctly assembled and closely matched.
- When assembling or disassembling, please be sure to assemble or disassemble each part according to the 'Assembly Method'.
- The materials used in this product meet the national safety standards, and it is normal to have residual odour.
- If you want to use the electric breast pump in the stage of less lactation, please use it under the guidance of doctors and nursing staff.
- Please don't over breast feed to avoid breast injury.
- If you feel pain caused by using this product, please stop using it immediately and ask the medical staff.
- Breast pump parts that come into direct contact with breast milk are personal care products, and mothers should not share them.
- This product is only for the purposes described in this manual.

WARNING

- Do not use it in the bath or shower.
- Be careful not to twist or knot the connecting tube.
- When not in use, please keep it out of the reach of children.

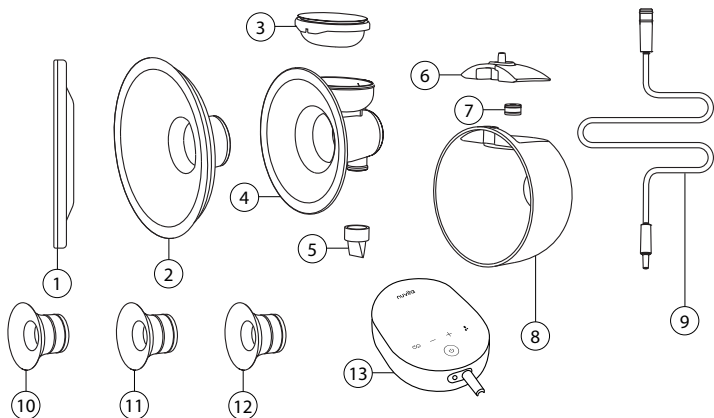
PACKAGING CONTENT

	Item	Quantity
1	Dust-proof cover	1
2	Breast cup	1
3	Suction cup	1
4	Breast pump body	1
5	Suction valve	1
6	Cap	1
7	Sealing ring for cap	1
8	Milk container	1
9	Connection tube	1
10	20 mm silicone reducer	1
11	22 mm silicone reducer	1
12	24 mm silicone reducer	1
13	Breast pump control unit	1
14	Velvet bag (not shown in the diagram)	1


ASSEMBLY METHOD

1. Insert the breast cup "A" to the breast pump body "D".
2. Put the suction bowl "C" on the top of the breast pump body. Ensure the edge of the suction silicone bowl is closely matched and the position limiting notch is installed in place.
3. Fix the suction valve "E" to the bottom of the breast pump body and ensure the suction valve is fully pushed in. The suction valve cannot be installed horizontally, and it must be installed as shown in the figure.
4. Fix the breast pump body assembly to the milk container "H" and make sure they are closely matched.
5. Push the cap assembly onto the milk container and make it fit together.
6. For extraction from a single breast, please connect the main breast pump unit and milk container assembly with a connection tube.
7. For double pumping (double breast pump mode), connect the main pump unit with a second set (Model 1287KDW Sold separately).


1

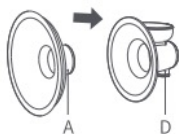
2

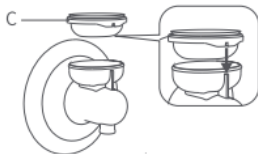3

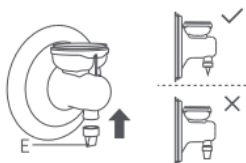
4


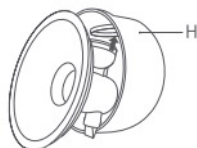
5

6

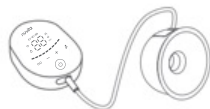7

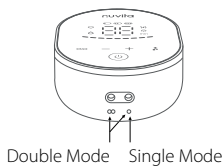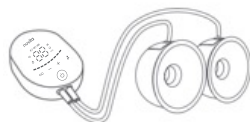Double Mode Single Mode

PRODUCT FEATURES

The Nuvita 1287CW Materno Smart Wear is designed to support breastfeeding mothers by mimicking a baby's natural nursing rhythm. This wearable electric breast pump is specifically engineered to help moms express breast milk efficiently and comfortably on the go, ensuring a seamless and pain-free pumping experience. It features 9 adjustable intensity levels, NipFit® functionality, and 3-phase Careflo® technology.

CAREFLO® TECHNOLOGY:

Unlike traditional 2-phase breast pumps, Careflo® introduces a third anti-stasis drainage phase alongside the standard stimulation and expression stages. This helps to:

- Promote a more consistent milk supply.
- Prevent clogged ducts and mastitis.
- Relieve breast tension and engorgement.
- Support a smoother, uninterrupted breastfeeding journey.

NIPFIT® FUNCTION:

Correction of flat or inverted nipples.

This function is designed to help correct flat or inverted nipples and support milk extraction through gentle, progressive stimulation:

- It helps the nipple protrude outward.
- It facilitates the baby's latch on to the breast.
- It supports the initiation and continuation of breastfeeding.

IMPORTANT: The NipFit® function only works in single mode.

CONTRAINDICATIONS

Do not use the NipFit® feature if:

1. You have open sores, cracks, lesions, or active infections on your nipple or areola.
2. You have had recent breast surgery.
3. You have a diagnosed condition of nipple inversion caused by disease or trauma, not congenital. If this is the case, consult your doctor first.
4. You experience sharp pain from the first application. Discontinue use immediately.
5. You have particularly sensitive skin or suffer from local dermatological conditions, such as eczema or psoriasis.
6. You are pregnant. Suctioning may stimulate contractions.

WARNING

The NipFit® feature is intended for postpartum use only. Do not use during pregnancy.

WARNINGS DURING USE


- Always start at the lowest level, Level 1. Do not increase the intensity until the nipple has gradually adapted, even if the sensation is comfortable.
- Do not exceed 5-10 minutes of continuous use of the NipFit® feature per session, especially in the first few weeks.
- Do not perform more than 2 correction sessions per day. Excessive use may cause irritation, swelling, or micro-tears to the nipple.
- Discontinue use immediately if you experience sharp pain, abnormal redness, persistent swelling, blisters, or discharge of fluid other than milk.
- This device is not a certified medical device for clinical therapy and does not replace consultation with a midwife, physician, or International Board Certified Lactation Consultant (IBCLC).
- A Grade 3, severely adherent inverted nipple requires medical evaluation before using any correction device.
- Prolonged, high-intensity use may cause swelling of the areola, temporarily making it more difficult for the baby to latch on. This can have the opposite effect to the desired one.

INSTRUCTIONS FOR SAFE USE:

- Make sure the silicone cup is the correct size, including the flange or reducers, before activating the NipFit® feature. A cup that is too large can amplify the suction on the areola in an uncontrolled manner.
- Before the session, apply a few drops of breast milk or pure lanolin cream to soften the area.
- During the NipFit® session, the nipple should be centered in the cup. An incorrect position can cause asymmetric suction and irritation.
- If you do not notice any improvement after 4-6 weeks of proper use, consult a healthcare professional.

BREAST PUMP CONTROL UNIT

- | | |
|---|--|
| A1. On and Off Switch button (long press to turn on) | B1. Level indicator |
| A2. Single/Double mode button | B2. Working time indicator(min) |
| A3. Working level decrease button | B3. Expression mode indicator |
| A4. Working level increase button | B4. Stimulation mode indicator |
| A5. Working mode switch button | B5. Drainage/Anti-stasis mode indicator |
| | B6. Single mode indicator |
| | B7. Minute indicator icon |
| | B8. Battery indicator |
| | B4+B5. NipFit® function activated |




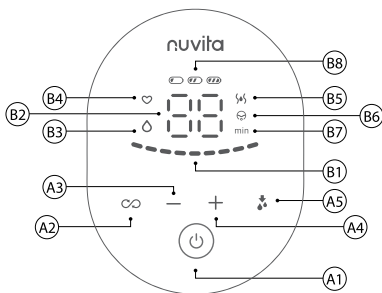
Intensity levels:

- | | | |
|---|-------|---|
|  | 1 - 9 | This is the Stimulation mode with fast and short suction rhythm. |
|  | 1 - 9 | This is the Expression mode. This simulates the suction of a baby. The rhythm of this phase is slower, deeper and longer. |
|  | 1 - 9 | This is the Drainage / Anti-stasis mode. This helps preventing Galactostasis. |

OPERATING INSTRUCTIONS

1. Wash hands before use, verify that the breast pump is perfectly clean and correctly installed according to the assembly instructions.
2. Align the center.
3. Adjust the nursing bra, fix the milk container, keep the body upright, and start expressing once the pump is properly connected to the main unit. (Note: do not lie down, bend, or move abruptly during extraction).
4. If the nipple diameter does not match the breast cup diameter, a silicone reducer can be used; insert the reducer into the breast cup and make sure it fits perfectly.
5. After finishing the extraction, remove the breast pump, pull out the cap from the milk container for disassembly, and pour the milk into the storage bottle or storage bag. (Note: do not exceed the maximum capacity of 180 ml).
6. After disassembly and cleaning, the parts can be placed in the sterilizer for disinfection. (Note: the tube and cap cannot be sterilized). Clean and sterilize the non-electronic parts for the next use.

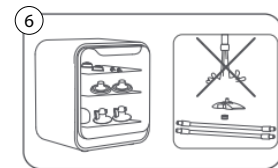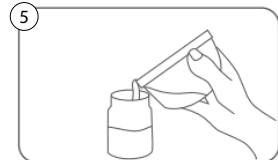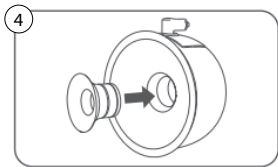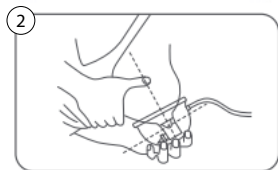BREAST PUMP OPERATION


- Power on, press **A1** power button. The screen of the machine will self test and enter stimulation mode automatically. The factory default level is level 1. Press **A4** to increase the level of suction intensity. After 2 mins of stimulation, the machine will auto switch to expression mode. Press **A4** again to increase the level of suction intensity. Machine will auto shut off after 30 mins working.
- Press **A5** to enter Expression mode if your milk starts to flow earlier. Press **A4** again to increase the level of suction intensity. Machine will auto shut off after 30 mins working.
- In Expression Mode, Press **A5** again to enter stimulation mode again at the same level of suction intensity in previous Stimulation session. After 2 mins of stimulation, the machine will auto switch to expression mode. Machine will auto shut off after 30 mins working.
- Press **A5** to enter Drainage mode when you prefer. Press **A4** again to increase the level of suction intensity. Machine will auto shut off after 5 mins.

INSTALLATION OF AC ADAPTOR

IMPORTANT: Please turn off the machine before connecting with the adaptor.


1. Insert the USB cable into the adapter port and plug it into the socket as shown in the figure. Insert the USB-C cable into the device.
2. When battery is low capacity, only the icon  will flicker, indicating that the machine needs to be charged.
3. When charging, the 3 battery icons on the screen will flicker.
4. When machine full charged, the icon  will be always on.

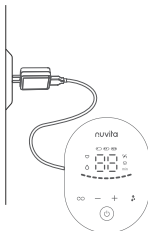HOW TO CHOOSE THE RIGHT REDUCER

The milk extraction should be a pleasant action and the size of the breast cup greatly affects the type of experience that is lived while using the breast pump. The Nuvita reducers are available in 3 sizes: 20, 22 or 24 mm.

1. With a ruler or a measuring tape, measure the diameter of the nipple (in mm), passing through the centre and excluding the areola.
2. Based on the mm measured, choose the breast reducer that is the closest to the measured value.

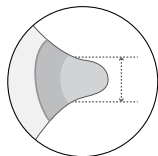**Size of
nipple**

UPTO
16 mm

UPTO
18 mm

UPTO
20 mm

**Size of
reducer**

Ø **20 mm**

Ø **22 mm**

Ø **24 mm**

CLEANING INSTRUCTIONS

- After each use, promptly clean the washable accessories.
- Do not use steel wool, sharp tools, or corrosive chemicals to clean this product.
- It is recommended to use soft sponge brushes and neutral detergent to clean this product.
- When disassembled and cleaned, accessories like the connector body, breast cup, milk container, suction cup, and valve can be placed in the sterilizer for disinfection. (Excluding the tube and cap group).
- The suction cup and valve are made of silicone. It is not advisable to sterilize them at high temperatures for long periods, as it may cause excessive aging or deformation and compromise their function.
- Do not immerse the connection tube or the cap in water. Clean them only with a clean, damp cloth. (If they get wet, they may cause permanent damage to the main unit).

- When disassembling and cleaning, pay attention to small components to avoid losing them, which may cause malfunction during assembly or operation.
 - While cleaning, gently wash the valve, suction cup, and other parts. If damaged, suction will be compromised and expression will not be possible.
 - This product cannot be heated or disinfected in a microwave oven.
-

WARRANTY - TERMS AND CONDITIONS

This product is covered by a 24-month warranty against manufacturing defects, starting from the date of purchase as evidenced by a receipt or invoice. If purchased with an invoice, the statutory warranty is reduced to 12 months. The 24-month warranty does not cover damage resulting from normal wear and tear, cosmetic damage, accessories or consumables (e.g. batteries, brush heads or parts subject to wear and tear). Batteries are covered by a 6-month warranty. Full warranty terms and conditions are available at [**www.nuvitababy.com/warranty**](http://www.nuvitababy.com/warranty)

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All rights reserved
Anteprima Brands International Ltd